

Ordering Information

Shape Type	Shape Code	Usable Length	Catalog No. 6Fr	Catalog No. 7Fr	Catalog No. 8Fr	
Judkins Left	JL3.5	100cm	HW60JL350P10000	HW70JL350P10000	HW80JL350P10000	
	JL3.5 ST	100cm	HW60JL350P10001	HW70JL350P10001	HW80JL350P10001	
	JL3.5 SH	100cm	HW60JL352P10000	HW70JL352P10000	HW80JL352P10000	
	JL3.5 STSH	100cm	HW60JL352P10001	HW70JL352P10001	HW80JL352P10001	
	JL4.0	100cm	HW60JL400P10000	HW70JL400P10000	HW80JL400P10000	
	JL4.0 ST	100cm	HW60JL400P10001	HW70JL400P10001	HW80JL400P10001	
	JL4.0 SH	100cm	HW60JL402P10000	HW70JL402P10000	HW80JL402P10000	
	JL4.0 STSH	100cm	HW60JL402P10001	HW70JL402P10001	HW80JL402P10001	
	JL4.5	100cm	HW60JL450P10000	HW70JL450P10000	HW80JL450P10000	
	JL4.5 ST	100cm	HW60JL450P10001	HW70JL450P10001	HW80JL450P10001	
	JL4.5 SH	100cm	HW60JL452P10000	HW70JL452P10000	HW80JL452P10000	
	JL4.5 STSH	100cm	HW60JL452P10001	HW70JL452P10001	HW80JL452P10001	
	Judkins Right	JR3.5	100cm	HW60JR350P10000	HW70JR350P10000	HW80JR350P10000
		JR3.5 SH	100cm	HW60JR352P10000	HW70JR352P10000	HW80JR352P10000
		JR4.0	100cm	HW60JR400P10000	HW70JR400P10000	HW80JR400P10000
JR4.0 SH		100cm	HW60JR402P10000	HW70JR402P10000	HW80JR402P10000	
Amplatz Left	AL0.75	100cm	HW60AL070P10000	HW70AL070P10000	HW80AL070P10000	
	AL0.75 SH	100cm	HW60AL072P10000	HW70AL072P10000	HW80AL072P10000	
	AL1.0	100cm	HW60AL100P10000	HW70AL100P10000	HW80AL100P10000	
	AL1.0 SH	100cm	HW60AL102P10000	HW70AL102P10000	HW80AL102P10000	
	AL1.5	100cm	HW60AL150P10000	HW70AL150P10000	HW80AL150P10000	
	AL1.5 SH	100cm	HW60AL152P10000	HW70AL152P10000	HW80AL152P10000	
	AL2.0	100cm	HW60AL200P10000	HW70AL200P10000	HW80AL200P10000	
	AL2.0 SH	100cm	HW60AL202P10000	HW70AL202P10000	HW80AL202P10000	
Short Amplatz Left	SAL0.75	100cm	HW60AL070P10002	HW70AL070P10002	HW80AL070P10002	
	SAL0.75 SH	100cm	HW60AL072P10002	HW70AL072P10002	HW80AL072P10002	
	SAL1.0	100cm	HW60AL100P10002	HW70AL100P10002	HW80AL100P10002	
	SAL1.0 SH	100cm	HW60AL102P10002	HW70AL102P10002	HW80AL102P10002	
	SAL1.5	100cm	HW60AL150P10002	HW70AL150P10002	HW80AL150P10002	
	SAL1.5 SH	100cm	HW60AL152P10002	HW70AL152P10002	HW80AL152P10002	
Amplatz Right	AR1.0	100cm	HW60AR100P10000	HW70AR100P10000	HW80AR100P10000	
	AR1.0 SH	100cm	HW60AR102P10000	HW70AR102P10000	HW80AR102P10000	
	AR2.0	100cm	HW60AR200P10000	HW70AR200P10000	HW80AR200P10000	
	AR2.0 SH	100cm	HW60AR202P10000	HW70AR202P10000	HW80AR202P10000	
Power Backup	PB3.0	100cm	HW60PB300P10000	HW70PB300P10000	HW80PB300P10000	
	PB3.0 SH	100cm	HW60PB302P10000	HW70PB302P10000	HW80PB302P10000	
	PB3.5	100cm	HW60PB350P10000	HW70PB350P10000	HW80PB350P10000	
	PB3.5 SH	100cm	HW60PB352P10000	HW70PB352P10000	HW80PB352P10000	
	PB4.0	100cm	HW60PB400P10000	HW70PB400P10000	HW80PB400P10000	
	PB4.0 SH	100cm	HW60PB402P10000	HW70PB402P10000	HW80PB402P10000	
Super Power Backup	SPB3.0	100cm	HW60SP300P10000	HW70SP300P10000	HW80SP300P10000	
	SPB3.0 SH	100cm	HW60SP302P10000	HW70SP302P10000	HW80SP302P10000	
	SPB3.5	100cm	HW60SP350P10000	HW70SP350P10000	HW80SP350P10000	
	SPB3.5 SH	100cm	HW60SP352P10000	HW70SP352P10000	HW80SP352P10000	
	SPB3.75	100cm	HW60SP370P10000	HW70SP370P10000	HW80SP370P10000	
	SPB3.75 SH	100cm	HW60SP372P10000	HW70SP372P10000	HW80SP372P10000	
	SPB4.0	100cm	HW60SP400P10000	HW70SP400P10000	HW80SP400P10000	
	SPB4.0 SH	100cm	HW60SP402P10000	HW70SP402P10000	HW80SP402P10000	
RCA Backup	RB0.75	100cm	HW60RB070P10000	HW70RB070P10000	HW80RB070P10000	
	RB0.75 SH	100cm	HW60RB072P10000	HW70RB072P10000	HW80RB072P10000	
	RB1.0	100cm	HW60RB100P10000	HW70RB100P10000	HW80RB100P10000	
	RB1.0 SH	100cm	HW60RB102P10000	HW70RB102P10000	HW80RB102P10000	
Special Curve	SC3.5	100cm	HW60SC350P10000	HW70SC350P10000	HW80SC350P10000	
	SC3.5 SH	100cm	HW60SC352P10000	HW70SC352P10000	HW80SC352P10000	
	SC4.0	100cm	HW60SC400P10000	HW70SC400P10000	HW80SC400P10000	
	SC4.0 SH	100cm	HW60SC402P10000	HW70SC402P10000	HW80SC402P10000	
Hockey Stick	HS01	100cm	HW60HS010P10000	HW70HS010P10000	HW80HS010P10000	
	HS01 SH	100cm	HW60HS012P10000	HW70HS012P10000	HW80HS012P10000	
	HS02	100cm	HW60HS020P10000	HW70HS020P10000	HW80HS020P10000	
Internal Mammary	HS02 SH	100cm	HW60HS022P10000	HW70HS022P10000	HW80HS022P10000	
	IM	100cm	HW60IM000P10000	HW70IM000P10000	HW80IM000P10000	
	IM SH	100cm	HW60IM002P10000	HW70IM002P10000	HW80IM002P10000	
Multi Purpose	MP01	100cm	HW60MP010P10000	HW70MP010P10000	HW80MP010P10000	
	MP01 SH	100cm	HW60MP012P10000	HW70MP012P10000	HW80MP012P10000	
	MP02	100cm	HW60MP020P10000	HW70MP020P10000	HW80MP020P10000	
	MP02 SH	100cm	HW60MP022P10000	HW70MP022P10000	HW80MP022P10000	



Number of side holes: 2
Not all shapes and sizes are available in all markets. Please check availability with your sales representative.

© 2016 ASAHI INTECC CO., LTD.
ASAHI, *ASAHI Hyperion* and *Corsair* are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.

Manufactured by:



ASAHI INTECC CO.,LTD.

3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan

Contact phone number: +81-561-86-9101

<http://www.asahi-intecc.com>

Distributed by:



Coronary Guide Catheter



Creating a stable environment for your PCI procedure



The guide catheter that can handle various sorts of difficult PCI conditions has arrived. The solution lies in here.

Inner Lumen

► The reasonable specifications for your stable PCI procedure

Enough room for your needs

- The large inner diameter diminishes the worry about friction between devices especially with 6F KBT, and provides more effective contrast visualization.



PTFE liner

- Smooth inner lumen reduces device insertion resistance.

	ID	OD
6 Fr:	1.80 mm (0.071")	2.09 mm
7 Fr:	2.05 mm (0.081")	2.40 mm
8 Fr:	2.28 mm (0.090")	2.70 mm

Flexible Tip

► Safe engagement and correct positioning made possible with our flexible-tip technology

Material



- Safe engagement is made possible due to the Urethane tip, the same material used for the tip of ASAHI Corsair, providing both flexibility and high visualization.



High visualization

All round processed tip

- Processed round tip walls result in atraumatic vessel engagement.

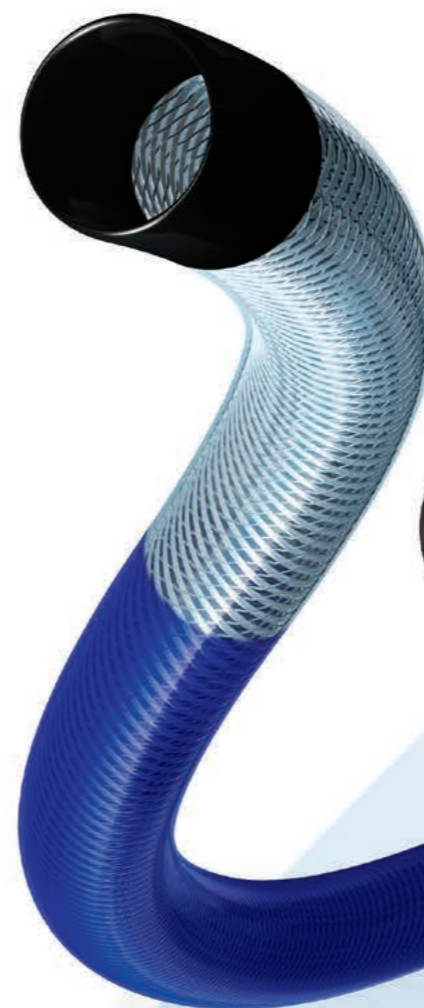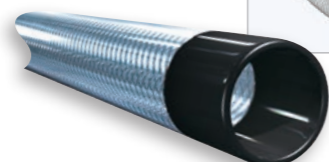

Backup

► Long-lasting backup even during challenging procedures

Hyper Shaft

- ASAHI INTECC's unique processing technologies prevent reduction of maneuverability and backup due to heat and moisture of blood during the procedure.


Reinforced 2nd curve

- A reinforced 2nd curve and excellent shape retention are key factors for high backup.

ASAHI INTECC-original 3rd curve

- ASAHI Hyperion's unique 3rd curve at the proximal part of the 2nd curve, improves backup by increasing the contact area between the catheter and the aortic wall.

*For JL, JLST, AL, SAL

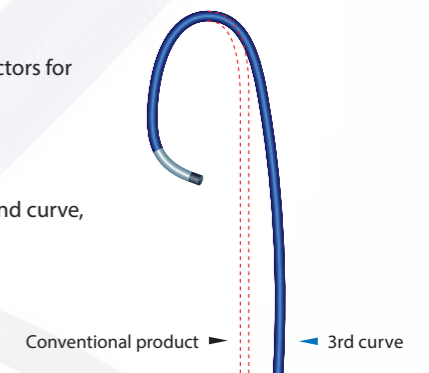


Maneuverability

► Provides superior maneuverability in tortuous vessels and prevents the catheter from jumping out of the ostium when engaging.

HENKA-Braid

- ASAHI INTECC's unique braiding technology provides well-balanced shaft with reinforced supportive shaft that gradually becomes softer towards the tip.

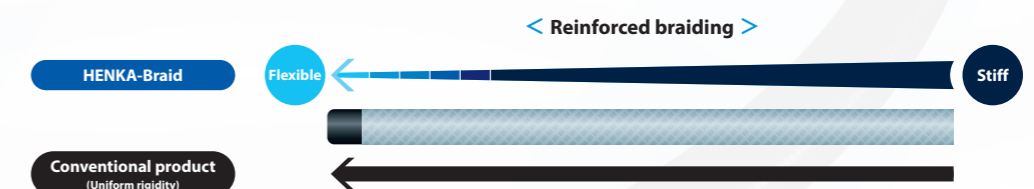

* HENKA is a Japanese word for transformation